

EXPERIMENT NO. 2

Aim: To determine the amount of Oxalic acid and H_2SO_4 present in 1 L. of solution **by** using N/10 NaOH and N/10 KMnO_4 solution.

Apparatus required: Burette, pipette, beakers, titration stand, conical flask and glass funnel etc.



Chemicals Required: NaOH, Oxalic acid, H_2SO_4 , Phenolphthalein and KMnO_4 .

Theory: This involves double titration.

1st Titration- Mixture (H_2SO_4 + Oxalic Acid) Vs NaOH:




NaOH reacts with oxalic acid as well as H_2SO_4 and find out the normality of oxalic acid and sulphuric acid.

Reactions involves

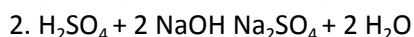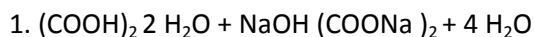


By titrating NaOH against the mixture, we find out the total normality of mixture. 2nd Titration-

Mixture (H_2SO_4 + Oxalic Acid) Vs KMnO_4 : The mixture of solution is treated with KMnO_4 solution which will react with oxalic acid in the presence of H_2SO_4 . By this titration we find out the normality and strength of oxalic acid.

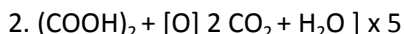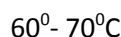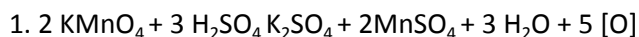

Complete reaction -

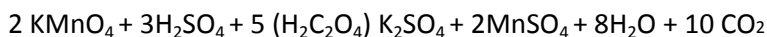

Procedure:

i) Titration of mixture (Oxalic Acid + H_2SO_4) and NaOH

1. Fill NaOH solution in the burette.
2. Now pipette out 10 ml of mixture (Oxalic Acid + H_2SO_4) in the conical flask.
3. Add 2-3 drops of Phenolphthalein indicator in the conical flask.
4. Now start pouring NaOH in the flask till the **colourless** solution in the flask turns to **pink**. And note the volume of NaOH used.
5. Repeat the experiment for concordant readings.

ii) Titration of Mixture (Oxalic Acid + H_2SO_4) and KMnO_4 :

1. Fill 50 ml of KMnO_4 solution in the burette. Since it is coloured solution we note upper meniscus

for taking the initial reading.

2. Now pipette out 10 ml of mixture (Oxalic Acid + H_2SO_4) in the conical flask and add 5 ml of H_2SO_4 in the flask.
3. Then heat the conical flask up to $60^\circ - 70^\circ\text{C}$.
4. Now start pouring KMnO_4 solution in the conical flask till the colourless solution in the flask turns to pink. And note the volume of KMnO_4 used.
5. Repeat the experiment for concordant readings.

Observations:

Titration i)

Solution in burette - NaOH

Solution in conical flask – (Oxalic Acid + H_2SO_4)

Indicator – Phenolphthalein

End Point – Colourless to Pink

Serial no.	Burette Reading		Volume Used(ml)
	Initial	Final	
1	0	28	28
2	0	29	29
3	0	29	29

Titration ii)

Solution in burette – KMnO_4

Solution in conical flask – (Oxalic Acid + H_2SO_4)

Indicator – KMnO_4

End Point – Colourless to Pink

Serial no.	Burette Reading		Volume Used(ml)
	Initial	Final	
1	0	10.2	10.2
2	10.2	20.4	10.2
3	20.4	30.6	10.2

Calculations:

Titration- i)

Mixture Vs NaOH

$$N_1 \times V_1 = N_2 \times V_2$$

(Mix.) (NaOH)

$$N_1 = N_2 \times V_2 / V_1$$

$$= 0.1 \times 29 / 10$$

$$N_1(\text{Mix.}) = 0.29 \text{ N}$$

Titration ii)

Mixture Vs KMnO₄

(Here, only Oxalic acid is titrated with KMnO₄)

$$N_1' \times V_1' = N_2' \times V_2'$$

(Oxalic acid) (KmnO₄)

$$N_1' = N_2' \times V_2' / V_1'$$

$$N_1' = 0.1 \times 10.2 / 10$$

$$N_1'(\text{Oxalic acid}) = 0.102 \text{ N}$$

$$\begin{aligned} \text{Normality of H}_2\text{SO}_4 &= N_1(\text{Mix.}) - N_1'(\text{Oxalic acid}) = \\ &0.29 - 0.102 \text{ N} \\ &= 0.188 \text{ N} \end{aligned}$$

$$\text{Strength of oxalic acid} = N_1' \times \text{Eq.wt.}$$

$$= 0.102 \times 63$$

$$= 6.426 \text{ g/L}$$

$$\text{Strength of H}_2\text{SO}_4 = 0.188 \times \text{Eq.wt.}$$

$$= 0.188 \times 49$$

$$= 9.212 \text{ g/L}$$

Result:


Oxalic Acid in the given mixture = 6.426g/L.

H₂SO₄ in the given mixture = 9.212 g/L


Structures:

Phenolphthalein

Colourless (Acidic Medium) Pink (Basic Medium)

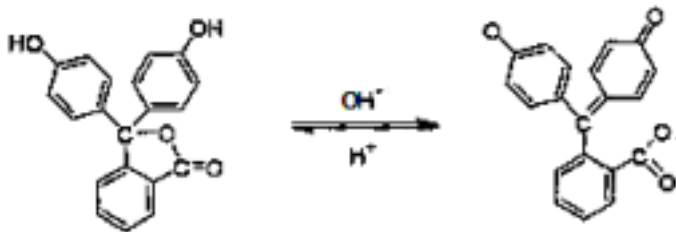

Structure of KMnO_4

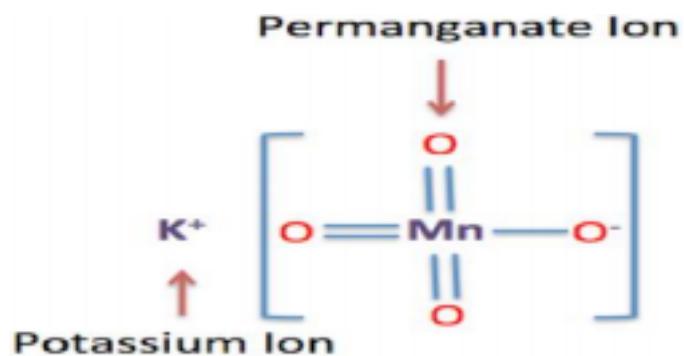

Precautions:

1. All the apparatus should be washed before use.
2. For measurement of coloured solution check the upper and lower meniscus of solution.
3. Rinse the burette with KMnO_4 solution and pipette with the given solution.
4. Always use freshly prepared KMnO_4 solution.

Applications:

1. It is used to determine the amount of various nutrients present in a particular food item.
2. It is used in innovation of medical sciences.