


EXPERIMENT NO. 6

Aim: To determine the alkalinity in given water sample.


Apparatus required: Burette, conical flask, pipette, beaker and measuring cylinder.

Chemical required: HCl, N/20 Na_2CO_3 solution, phenolphthalein, methyl orange.

Theory: Alkalinity is a measure of the acid buffering capacity of water. Alkalinity of a sample of water is due to the presence of hydroxide ion (OH^-) bicarbonate ion (HCO_3^-) and carbonate ion (CO_3^{2-}) or the mixture of any of the two ions present in water. But the (OH^-) and (HCO_3^-) ions together is not possible as they combine together to form (CO_3^{2-}) ion.

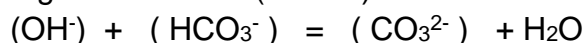

Reactions involved:

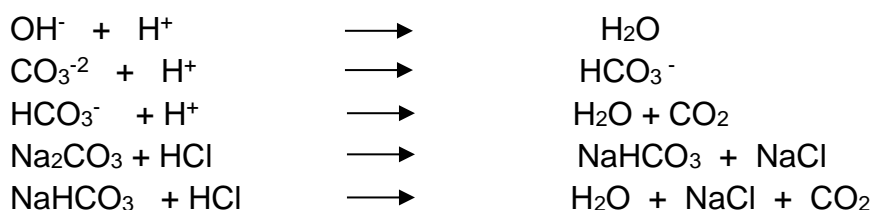

Procedure:

I- Standardization of HCl

1. 50 ml HCl is filled in the burette.
2. 10 ml Na_2CO_3 solution is pipette out in the conical flask.
3. 2-3 drops of Methyl Orange is added and titration done is till end point observed.

II- Titration with H_2O Vs. HCl

1. 50 ml HCl is filled in the burette.
2. 10 ml H_2O is pipette out in the conical flask.
3. 2-3 drops of Phenolphthalein is added.
4. Run the HCl in conical flask.
5. Titration is done till the solution turns pink to colourless.
6. 2-3 drops of Methyl Orange is added and titration is continued till the Orange-Pink colour is observed.

Observations:

Standardization of HCl

Solution in burette = HCl
Solution in conical flask = 10 ml of Na_2CO_3
Indicator = Methyl Orange
End Point = Pink colour

Observation Table:

Serial no.	Burette Reading		Volume of HCl (ml)
	Initial	Final	

1	0	17.9	17.9
2	17.9	35.8	17.9

HCl Vs H₂O

Solution in burette = HCl

In conical flask = Water (10 ml)

Indicator = Phenolphthalein and Methyl Orange

End Point 1st = Pink to colourless

End Point 2nd = Colourless to pink

Observation Table:

Serial no.	Phenolphthalein			Methyl orange			V ₁ + V ₂ (ml)
	Initial	Final	V ₁ (ml)	Initial	Final	V ₂ (ml)	
1	0	16.5	16.5	16.5	20.5	4.0	20.5
2	20.5	37	16.5	37	41	4.0	20.5

Calculations :

i)Standardization of HCl

$$\begin{aligned}
 N_1 V_1 &= N_2 V_2 \\
 (\text{HCl}) & \quad (\text{Na}_2\text{CO}_3) \\
 N(\text{HCl}) &= (1/20) \times 10 / 17.9 \\
 N(\text{HCl}) &= 0.027 \text{ N}
 \end{aligned}$$

II)HCL Vs Water

a) For Phenolphthalein-

$$\begin{aligned}
 N_1 V_1 &= N_2 V_2 \\
 (\text{H}_2\text{O}) &= (\text{HCl}) \\
 N(\text{H}_2\text{O}) &= (0.027 \times 16.5) / 10 \\
 &= 0.033 \text{ N} \\
 \text{Strength} &= 0.033 \times 50 \\
 P &= 1.65 \text{ g/L} \\
 &= 1650 \text{ ppm}
 \end{aligned}$$

b) For Methyl Orange-

$$\begin{aligned}
 N_1 V_1 &= N_2 V_2 \\
 (\text{H}_2\text{O}) &= (\text{HCl}) \\
 N(\text{H}_2\text{O}) &= (0.027 \times 20.5) / 10 \\
 &= 0.055 \text{ N} \\
 \text{Strength} &= 0.055 \times 50 \\
 M &= 2.76 \text{ g/L} \\
 &= 2767 \text{ ppm}
 \end{aligned}$$

$P > (1/2) M$, so OH^- and CO_3^{2-} ions are present in water


Concentration of $\text{OH}^- = (2 P - M) = 533 \text{ ppm}$

Concentration of $\text{CO}_3^{2-} = 2 (M - P) = 2234 \text{ ppm}$

Structure:


Methyl Orange

Pink (Acidic Medium)

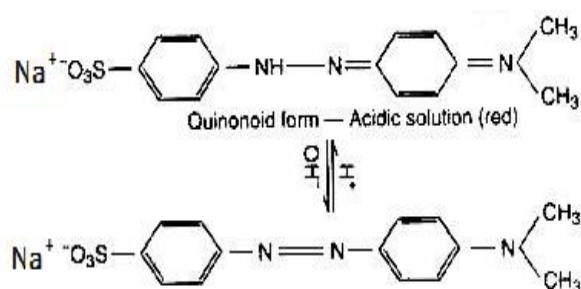

Yellow (basic Medium)

Phenolphthalein

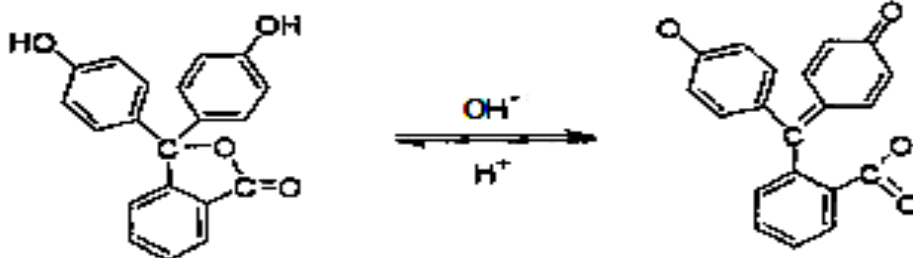

Colourless (Acidic Medium)

Pink (Basic Medium)

Result: Since $P > (M/2)$, thus OH^- and CO_3^{2-} ions are present in the sample of water.

OH^- Concentration = 1250 ppm

CO_3^{2-} Concentration = 800 ppm

Precautions:

1. Wash the apparatus before use.
2. Rinse the burette with the solution to be filled in it.
3. For measuring coloured solution the upper meniscus is taken and for colourless solution measure the lower meniscus is taken.